REMINGTON®

i-LIGHT® REVEAL™

HAIR REMOVAL

INTENSE PULSED LIGHT (IPL) HAIR REMOVAL
FOR THE FACE AND BODY

GUARANTEE ⁺¹
2 yr
IPL2000

SKIN CHART

THE SKIN CHART

THE HAIR CHART

WARNING : i-Light REVEAL is not suitable for very dark or black skin, white/grey or very light blonde/red hair.

For facial use: Only suitable for female facial use below the cheekbone.

DO NOT use on or around the eyes, eyebrows or eyelashes. Doing so can cause serious and permanent eye injury.

Thank you for buying your new Remington® product.

Please read these instructions carefully and keep them safe. Remove all packaging before use.

IMPORTANT SAFEGUARDS

- 1 **WARNING – TO REDUCE THE RISK OF BURNS, ELECTROCUTION, FIRE, OR INJURY TO PERSONS:**
- 2 Use this appliance only for its intended use described in this manual.
- 3 Do not use if it is damaged or malfunctions.
- 4 Do not leave the appliance unattended while plugged in.
- 5 Keep the power plug and cord away from heated surfaces.
- 6 Do not use the product with a damaged cord. A replacement can be obtained via our International Service Centres.
- 7 Use and store the product at a temperature between 15° and 35°C.
- 8 Only use the parts supplied with the appliance.
- 9 Do not twist or kink the cable and don't wrap it around the appliance.
- 10 The body of this appliance is not washable or water resistant. Do not put the appliance in liquid; do not use it near water in a bath tub, basin or other vessel; and don't use it outdoors.
- 11 Make sure the power cable and adaptor do not get wet. Do not reach for your unit if it has fallen into water - immediately unplug it from the electrical outlet.
- 12 Do not plug or unplug the appliance with wet hands.
- 13 This appliance can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- 14 Keep the appliance and cable out of reach of children.
- 15 This appliance should be supplied with approved safety isolating adaptors SW-120060EU (for Europe) and SW-120060BS (for United Kingdom) with the output capacity of 12 dc; 600mA (adaptor output).
- 16 DO NOT cover any part of the ventilation slots during use, as they help to keep the device cool during operation. This could result in the device overheating.
- 17 DO NOT use if the box packaging or the device appears damaged or if any part of the device is cracked or broken. Doing so could increase the risk of skin or eye injury from the flash lamp.
- 18 This appliance is not intended for commercial or salon use.

i-LIGHT REVEAL KEY FEATURES

 KEY FEATURES**DEVICE**

- 1 Flash Window
- 2 Treatment Head
- 3 Treatment Head Release Tabs
- 4 Skin Tone Sensor
- 5 Flash Button
- 6 LCD Display
- 7 Lint Free Cloth
- 8 Storage Pouch (not shown)

LCD DISPLAY

- 9 Battery Life
- 10 Session Flash Counter
- 11 Bulb Life
- 12 Energy Level
- 13 Energy Level Button
- 14 Power On/Off
- 15 Power Adaptor

👁 GETTING ACQUAINTED WITH i-LIGHT REVEAL

★ What is i-Light?

The i-Light system is a revolutionary light-based device designed for the removal of unwanted body hair and female facial hair at home.

i-Light REVEAL utilises similar clinically-proven technology that professional dermatologists and salons use, optimised for personal treatment of smaller areas.

★ What is Intense Pulse Light (IPL) and how does the i-Light REVEAL work?


IPL - (Intense Pulsed Light) is a well-established technology used worldwide for over 15 years. i-Light REVEAL uses this technology to flash your skin with an intense pulse of light.

The pulse works under the skin's surface to target hairs at the root, helping to stop their growth. (Fig 1)

The energy in the pulse of light targets dark colouring, and best results are seen on darker hair.

To avoid targeting (and possibly damaging) dark skin, you should only use i-Light REVEAL on light to medium skin tones (I-IV).

Fig. 1

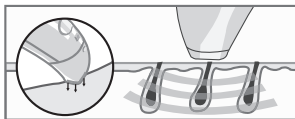★ What to expect from i-Light REVEAL

Results aren't immediate, but most users will begin to notice fewer hairs within 1-2 weeks after beginning treatment. During this time it may look like the hair is still growing, but it is actually being 'pushed out'.

Fig. 2.

Fig. 2

It's important to continue with the full treatment regime to ensure that every active hair is treated. A full treatment regime is every 2 weeks for three treatments.

Individual results may vary, multiple treatments may be required for optimal results.

Some hair growth will occur that is less dense, finer and lighter in colour than the original hair. This regrowth is perfectly normal and should be expected.

After completing the initial treatment regime, continuous regular weekly use will give good results within 6-12 weeks (darker skin tones may take longer).

WARNINGS AND SAFETY PRECAUTIONS

◆ IMPORTANT SAFETY CAUTIONS

Before you start using the i-Light REVEAL:

Be sure to read all Warnings and Safety Information. Before you begin, check to see if i-Light REVEAL is suitable for you.

Use the skin chart provided on the box and at the front of this manual and the integrated skin sensor to determine if this device is right for you.

Consult the hair chart at the front of this manual to ensure the device is right for you.

★ Skin type

See the skin colour chart on page 2.

- Do not use on naturally dark skin, as it may result in burns, blisters and changes in skin colour.
- Do not use on tanned skin or after recent sun exposure, as it may cause burns or skin injury.
- i-Light REVEAL is not effective on naturally white, grey, blonde or red body or facial hair.

★ Areas not to treat:

- Do not use on the scalp or ears.
- Do not use on male face or neck.
- Do not use above the cheekbone (female facial use).
- Do not use on nipples, areola or genitals.
- Do not use if you have tattoos or permanent makeup in the area to be treated.
- Do not use on dark brown or black spots such as moles, birth marks or freckles.
- Do not use on an area of recent surgery, deep peel, laser resurfacing, scars or skin that has been damaged with burns or scalds.

★ When not to use/When to avoid using i-Light REVEAL

- DO NOT use if you have TANNED SKIN in the intended treatment area, as your skin may be too dark and treatment in tanned areas could temporarily cause skin swelling, blisters, or could change the colour of your skin. Wait until your tan has faded before treating.
- DO NOT use on skin with cosmetics applied, as treating in those areas could temporarily cause skin swelling, blisters or could change the colour of your skin. Make sure to thoroughly clean the intended treatment area prior to use.
- DO NOT flash the same treatment spot multiple times during the same session, as this could temporarily cause skin swelling, blisters, or could change the colour of your skin.
- DO NOT treat any area more than once per week, as this will not speed up results, but could increase likelihood of skin swelling, blistering or changing the colour of your skin.
- DO NOT use on damaged skin, including open wounds, cold sores, sunburns, or areas where you have had a skin peel or other skin resurfacing procedure within the last 6-8 weeks.
- DO NOT use on irritated skin, including rashes or swollen skin. These conditions can make the skin more sensitive. Wait for the affected area to heal before use.

- DO NOT use if your skin is sensitive to light. In order to determine if you can use i-Light REVEAL, perform a test on a small patch of skin and wait 24 hours to determine if there are any adverse reactions before treating more broadly.
 - DO NOT use if you are or may be pregnant or are breast-feeding. i-Light REVEAL has not been tested on pregnant women. Hormonal changes are common during pregnancy or breast-feeding and could increase skin sensitivity and the risk of skin injury.
 - DO NOT flash against any material other than skin.
 - DO NOT use this model (IPL2000) for any purpose other than hair removal.
 - DO NOT treat any area where you may someday want your hair to grow back. The results may be permanent and irreversible.
 - ★ CAUTION: Avoid treating areas that contain only vellus hairs (very thin, fine hair known as “peach fuzz”) as this may result in undesired hair re-growth.
 - DO NOT use on the face if you have a history of cold sores; use of IPL or any treatment that irritates skin can cause a recurrence of cold sores.
 - If you’ve had professional or clinical or self-administered laser or IPL or electrolysis hair removal treatments, wait at least 6 months before using i-Light REVEAL. If you’ve waxed, plucked or epilated hair in the treatment area, wait at least 6 weeks before using i-Light REVEAL. The device uses the hair under your skin to disable your hair at the root, so it is not effective if the hair has been pulled out.
 - DO NOT use directly over or near Active Implantable Medical Devices such as implantable pacemaker, internal defibrillator, or insulin devices. Although no problems associated with the use of this device have been reported during clinical studies, it is recommended to observe all active implantable medical device warnings regarding the use of electronic devices to prevent any failure or interruption of an active implantable device.
 - Do not use over or near anything artificial such as silicon implants, Implanon contraceptive implants, pacemakers, subcutaneous injection ports (insulin dispenser) or piercings.
 - Wireless communications equipment may affect the i-Light REVEAL device and should be kept at least 3.3 m away.
- Do not use over or near metal filings, dental implants or braces.

PREPARING FOR USE

★ INSTRUCTIONS FOR CHARGING

- Before your first use; ensure the device is fully charged.
- Plug the power adaptor into the device and then into the socket.
- One full charge will take up to 1.5 hours.
- Dependent upon the energy level selected; one full charge will provide approximately 100 flashes (Lower energy levels will provide more flashes per charge).

★ Quick Charge Function

- The device is also equipped with a quick charge function.
- To utilise this function, charge the device for 15 minutes.
- This will provide approximately 20 flashes.

- When the battery is charging, the LCD Display will show the charging icon Fig. 3
- When the battery is fully charged, the LCD Display will show the battery full icon Fig. 4

Fig. 3

Fig. 4

- During use the status of the battery is shown on the LCD display by the 'Battery Life' icon. Fig. 5


Fig. 5

Battery Life

Low Battery Life

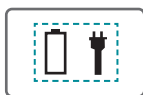Battery Life Empty

- * NOTE: The device cannot be operated whilst it is plugged in or being charged.
- * TIP: We recommend charging the device after each use to ensure that the device is fully-charged and ready for your next treatment.

👁️ USING i-LIGHT REVEAL FOR TREATMENT

* Skin Tone Test

1. Review the Warnings and Safety Precautions.
2. Familiarise yourself with the features of the i-Light.
3. Consult the Skin Tone Chart to ensure your skin colour is in the acceptable range. Fig. 6




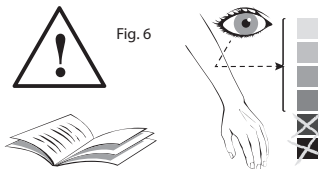
Fig. 6

- Turn on the device by pressing the power button. Fig. 7

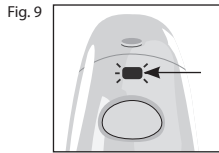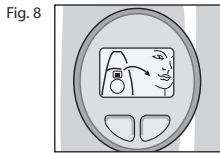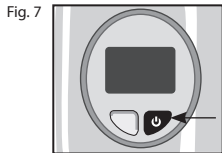- NOTE: Each time the device is turned on, you will be prompted to test your skin tone. The LCD display will be showing the 'Skin Tone Sensor' icon. Fig 8
 - The Skin Tone Sensor is located on the back of the device and will show a flashing orange light. Fig. 9
- Firmly press the Skin Tone Sensor (flashing orange light) onto the skin of the intended treatment area and hold for 2 seconds. Fig. 10.

Fig. 10


- The reading will take place automatically.
 - Solid orange light and 'buzz' = Skin may be too dark
 - Solid green light and 'beep' = Skin tone is OK
- NOTE: You must repeat the skin tone test for each area you wish to treat.

Energy Level Selection

- The device will automatically be set to energy Level 1 each time the device is powered ON. To change the level, press the Energy Level Selection Button. Fig. 11

Fig. 11

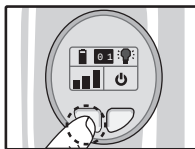- After reaching level 3, pressing the Energy Level Button again will reset back to energy level 1.
 - The i-Light REVEAL is equipped with 3 energy levels. Level 1 (1 bar) is the lowest and level 3 (3 bars) is the highest.
 - ✳ **TIP:** For the most effective results, always use the highest energy level that does not cause discomfort on the skin.
 - ✳ **Test the i-Light REVEAL on your skin**
 - Before your first full treatment, test the i-Light REVEAL device on a small patch of skin and wait 24 hours to ensure there are no adverse reactions.
 - Start at energy level 1 and flash once.
- Flash the skin by placing the flash window firmly against your skin. You will hear a 'beep beep' and the skin tone sensor should show a solid green light. Press the flash button to deliver a flash. Fig. 12

Fig. 12

- ✳ **TIP:** It is helpful to use a mirror to see the skin tone sensor. If the reflected light is unpleasant, you may wish to close your eyes before each flash.
- Mild to moderate discomfort is normal. If a flash delivered at Energy Level 1 feels tolerable, increase to Energy Level 2. Move the device to a different spot and flash the i-Light REVEAL one more time. Repeat this process until you find the setting that feels right for you.
- Wait 24 hours and then look at the test area. If the skin appears normal, proceed with your first full treatment at that the energy level you determined is right for you.
- ✳ **NOTE:** You must test the device on each area of the face or body you wish to treat.
- ✳ **Treat desired area(s) with the i-Light device**
- ✳ **Prepare your skin for treatment**
 1. Ensure the area to be treated is clean and free from oils, deodorants, perfume, make up, lotions and creams. Shave or trim the hair from the area to be treated, and ensure the skin is clean and dry. If trimming; the hair must be trimmed to less than 0.5mm. Fig. 13

Fig. 13

- Never use wax, epilation, tweezers or depilation products to remove the hair because they counteract the IPL process.

✱ Prepare the device for treatment

Use the skin tone sensor to unlock the unit. Fig. 14
Select Energy Level Fig. 15

Fig. 14


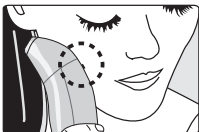
Fig. 15

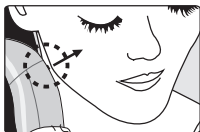👁 TREATMENT


- Place the flash window firmly against the skin. The light will be solid green if the flash window is touching the skin. Fig. 16

Fig. 16


- Press the Flash Button to deliver a flash. You should expect to see a bright flash of light (similar to a camera flash), and you should feel a mild sensation of warmth.
- ✱ NOTE: The device can flash approximately every 5 seconds. During this time an hour-glass icon will appear on the display Fig. 17. When the device is ready to flash again; the icon will disappear and the light will be solid green (if treatment head is in full contact with skin).

Fig. 17

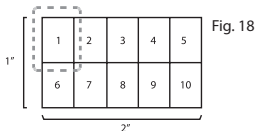
- Treat the skin in a grid-like pattern. Fig. 18

- For best results, please follow the guidelines in fig.19 which show the recommended number of flashes per facial area.

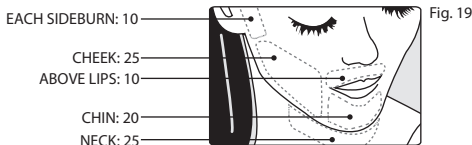- TIP:** Stretch Skin for Better Skin Contact – the contours of the skin may make it difficult for the skin contact sensors to be fully depressed in certain areas. You may find it easier to try and wrap your upper lip over your teeth or use your tongue to stretch your cheeks for better skin contact.
- To help you keep count of the number of flashes you have used in one treatment; the device is equipped with a 'Session Flash Counter'. Fig. 20. This will reset each time the device is switched on.

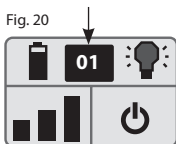POST-TREATMENT CARE

After treatment, you may experience slight redness or a warm sensation on your skin. This is normal and will disappear quickly. To avoid irritation to your skin after a treatment, take the following precautions:

- Avoid sun exposure for 24 hours after a treatment. Protect the skin with SPF 30 for 2 weeks after each treatment.
- Do not prolong sun exposure such as sunbathing, using a tanning bed, or self-tan for at least 2 weeks after the last treatment.
- After treatment, keep the area clean and dry and drink plenty of water to keep skin hydrated.
- Do not take hot baths, showers, or use steam rooms and saunas for 24 hours after treatment.
- Do not take part in contact sports for 24 hours after treatment.
- Do not depilate (waxing, plucking, threading or creams) during the treatment – shaving is acceptable as long as you avoid shaving 24 hours after each treatment.
- Do not use bleaching creams or perfumed products for 24 hours after treatment.


CLEANING YOUR i-LIGHT REVEAL DEVICE

- ◆ **CAUTION:** Before cleaning your i-Light REVEAL, make sure that the device is switched off and the charging cord is disconnected from the device.
- Regular cleaning helps to ensure optimal results and a long life for the device.
- Use a clean, lint-free cloth (one is provided with the device) or a cotton swab to gently wipe the surface and the Flash Window, Treatment Head and the Skin Tone Sensor. If needed, you can dampen the cloth or cotton swab with water and mild soap.
- Allow the device to dry completely before using again.
- ◆ **WARNING:** If the Flash Window is cracked or broken, the unit must not be used. Never scratch the filter glass or the metallic surface inside the Treatment Head.
- ◆ **CAUTION:** The i-Light is a high voltage device. Never immerse in water. Never clean the unit or any of its parts under a tap or in a dishwasher.
- Do not use petroleum-based or flammable cleaning agents because of the risk of fire. Never use scouring pads, abrasive cleaning agents or aggressive liquids such as oil or acetone to clean the unit.

TROUBLESHOOTING / STORAGE / MAINTENANCE

i-Light Reveal Device Maintenance

- ◆ **CAUTION:** Before performing maintenance on your device, ensure the device is switched off and the power cord is disconnected from the base unit.
- ◆ **Replacing the Treatment Head**
 - The treatment head can deliver 6,000 flashes regardless of the energy level setting. When the Treatment Head has used up its flashes, it must be replaced.
 - The 'Bulb Life Remaining' icon on the display warns you when the number of remaining flashes is low. Fig. 21.

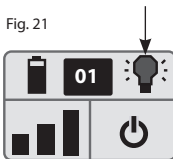- Bulb Life Remaining Display Fig. 22.

Fig. 22

Status	Full	Partially Full	Empty
Bulb life remaining display	A new bulb provides 6,000 flashes and the Bulb Fill Colour will be solid	When there are fewer than 300 flashes left, the Bulb Fill Colour will turn to black	Bulb Life is empty and the entire bulb is flashing. The device will 'beep'. Treatment Head must be replaced
ICON			

✱ NOTE: The Treatment Head must be replaced if the Lamp Tube is broken.

1. Press the treatment Head Release Tabs and gently pull to remove the Treatment Head. Fig. 23.
2. Replace with a new Treatment Head, ensuring it snaps into place. Fig. 24.

Fig. 23


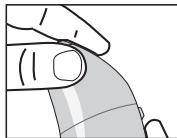
Fig. 24

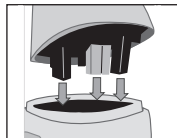✱ Storage

- Switch off the unit, unplug it and let it cool down for 10 minutes before storage.
- Store the unit in a dry place at a temperature between 15° C and 35° C.

✱ Troubleshooting

- Always read these instructions fully before using i-Light REVEAL.
- Refer to this troubleshooting guide if you experience any problems with i-Light REVEAL, as this section addresses the most common problems you could encounter.
- If you have followed the instructions in this section and continue to experience problems, please contact your local Remington® Service Centre for further assistance.

Problem	Possible Cause	Action
The device will not switch on.	The device isn't charged. You should hear the device 'beep' twice and the LCD display will show this screen. 	Ensure the device is fully charged - plug the charging cord into the wall socket. The Battery Charging icon should appear on the LCD display. When the device is fully charged, the Battery Life indicator will be full: 
	The Treatment head is not fully attached. You should hear two 'beeps' and the LCD display will show this screen. 	Ensure the Treatment Head is properly attached to the device. When the Treatment Head is attached correctly, you should hear a 'click' and the lights and 'beeping' will stop.

Problem	Possible Cause	Action
When I hold the device against my skin, nothing happens. The device doesn't 'beep' even after a few seconds.	You may need to perform a Skin Tone Test. The LCD display will show this screen: 	Place the flashing orange light against the skin of the intended treatment area. The device may require another Skin Tone Test if: -The device is not used for five minutes - You take a break mid-treatment - The device is switched off and then on again
	The Treatment Head is not in full contact with the skin.	Check that the Treatment Head is touching your skin. The device will not flash unless the Treatment Head is making full contact with the skin.
	You are trying to treat a small area such as between the eyebrows or near the eyes.	For your safety, the device should NOT be used anywhere near the eyes, eyebrows or eyelashes.

Problem	Possible Cause	Action
A flash of light is not emitted when I press the Flash Button	The Treatment Head is not in full contact with the skin.	Check that the Treatment Head is touching your skin. The device will not flash unless the Treatment Head is making full contact with the skin.
	You have not waited enough time between flashes.	The device requires approximately 5 seconds to recharge between flashes. When the device is ready to flash again, the light will become solid green.
	The Treatment Head needs to be replaced; the LCD display will show the 'Bulb Life Empty' icon and will be flashing on and off.	Ensure the Bulb Life Icon is not empty and flashing. If it is flashing; turn the device off and replace the Treatment Head.

Problem	Possible Cause	Action
A flash of light is not emitted when I press the Flash Button	The device needs some time to cool down. You should hear the device 'beep' once and the LCD display will show this screen: 	Wait a few minutes for the device to cool down. Do not turn the device off as this screen will disappear once the device is ready to use again. Under normal conditions, this is uncommon, however this can arise if the unit becomes too hot.
	The room is too warm for the device to cool down.	Use i-Light REVEAL in a cooler environment where the temperature is less than 27°C (or 80°F).

Problem	Possible Cause	Action
I performed a Skin Tone Test and the light turned orange and 'buzzed'.	Your skin may be too dark for treatment. Refer to the Skin Chart at the beginning of this manual.	You cannot use the i-Light REVEAL device with your skin tone. Please contact your local service centre.
	You didn't use the Skin Tone Sensor correctly.	If the Skin Tone Sensor (flashing orange light) isn't firmly touching the skin, the device will make a rejection tone. Try pressing the Skin Tone Sensor firmly against the skin.
	The Skin Tone Sensor isn't working.	Please contact your local service centre.

Problem	Possible Cause	Action
Dirt or debris is stuck to the Flash Window.	The device needs to be cleaned.	Follow the cleaning instructions provided in this manual.

Problem	Possible Cause	Action
The device: - Isn't working - Is beeping - Is showing this symbol (including a reference number): 	The device isn't working properly and may need to be serviced.	Contact your local service centre; and provide the reference number from the display.

- ⚡ **DO NOT** attempt to open or repair your device. Only authorised i-Light REVEAL repair centres are permitted to perform repairs.
- Opening the device may expose you to dangerous electrical components and to pulsed light energy, either of which may cause serious bodily damage and/or permanent eye injury. Trying to open the device may also damage the device and will void your guarantee. Please contact your local service centre if you have a device in need of repair.

ENVIRONMENTAL PROTECTION

To avoid environmental and health problems due to hazardous substances in electrical and electronic goods, appliances marked with this symbol must not be disposed of with unsorted municipal waste, but recovered, reused or recycled.

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© 2014 SBI
Remington® is a Registered Trade Mark of Spectrum Brands, Inc.,
ÜRETİCİ FİRMA BİLGİSİ: VARTA Consumer Batteries GmbH & Co. KGaA, Alfred-Krupp-Str. 9, 73479 Ellwangen, Germany.
VARTA Consumer Batteries GmbH & Co. KGaA, Alfred-Krupp-Str. 9, 73479 Ellwangen,